

# Lightweight wheelchair Freetec

Item No. 900500500



- Simple adjustment from short to long armrests
- Front fork camber adjustment
- Footrests adjustable in five steps from 0 to 32 degrees
- Wheelbase extension possible
- Available in eight seat widths
- Handle height adjustable to four positions
- Back cover depth and height adjustable
- Armrest reinforced with internal steel insert for greater stability
- Heel straps provide additional hold
- 200 x 50 mm PU steering wheels and 24" x 1 3/8" full-rubber drive wheels
- Extra-reinforced front wheel forks
- Three-position seat height adjustment from 45 cm to 50 cm in the standard model for maximum comfort
- Continuous adjustment of leg length
- Quick-release axles as standard for fast wheel installation or removal
- Optionally available with drum brake for escort
- Crash test according to ISO 7176-19: approved for passenger transport in vehicles in combination with all approved passenger and wheelchair restraint systems
- The Freetec with common part principle: gain flexibility and reduce your stock level – with the common part principle: many Freetec replacement parts and accessories are compatible with all Litec 2G, Litec 2G Plus, Rotec, Ecotec 2G, and Multitec wheelchairs

## Technical Data

# Lightweight wheelchair Freetec

Item No. 900500500



|                 |              |                        |                |
|-----------------|--------------|------------------------|----------------|
| Order number    | 900500500    | Armrest height         | 21.5 cm        |
| Colour          | black-gray   | Width between armrest  | 0 cm           |
| Drum brake      | No           | Footrest adjustment    | 34.5 - 48.5 cm |
| Overall width   | 66 cm        | Drive wheel dimensions | 24" x 1 3/8"   |
| Overall length  | 103 cm       | Front wheel dimensions | 200 x 5 cm     |
| Overall height  | 92 - 99.5 cm | Folded width           | 28 cm          |
| Seat width      | 48 cm        | Total weight           | 16.6 kg        |
| Seat depth      | 43 - 50 cm   | Max. Load              | 125 kg         |
| Seat height     | 45 - 50 cm   | PU Pal.                | 8 Pal.         |
| Backrest height | 38.5 - 46 cm |                        |                |

# Lightweight wheelchair Freetec

Item No. 900500500

**drive**

**DeVilbiss**  
HEALTHCARE

## Accessories



**Power knot system Prosafe**

Item No.900100783




**Height- and length-adjustable armrest**

Item No.900300100




**Brake lever extension**

Item No.900300110

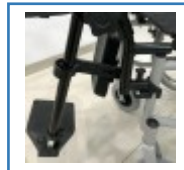

**Folding and lockable leg supports**

Item No.900300170

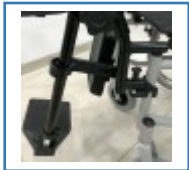


**Folding and lockable leg supports**

Item No.900300171



**Folding and lockable leg supports**

Item No.900300172

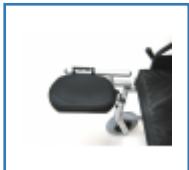

**Stump support**

Item No.900300180



**Stump support**

Item No.900300190



**Cane holder for wheelchairs**

Item No.900300311



**Lap belt L**

Item No.900300320



**PU wheel 150 x 35 mm with wheel fork**

Item No.900300340



**Oxygen bottle holder for wheelchairs**

Item No.900300350



**Wheelbase extension adapter block  
Freetec**

Item No.900300360



**Wheelchair one-hand drive for Rotec  
and Freetec**

Item No.900100778

# Lightweight wheelchair Freetec

Item No. 900500500

**drive**

**DeVilbiss**  
HEALTHCARE



**Seat cushion**

Item No.900300272